


*Anytime, anywhere, anyone.*

## FEATURES



**PACE-PULSE DETECTION AND REJECTION FUNCTION.** The A16 filters external pulse, like a pace-maker, and removes from the ECG reading.

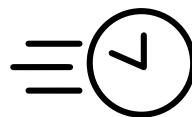

**FAST SHOCK DELIVERY**  
Shock delivery within 10 seconds.



**MULTILINGUAL**  
English, Chinese (Mandarin), Korean.



### OPERATION

**Fully-automatic AED** - no need to press the shock button. **Semi-automatic AED** - flashing button indicates 'ready for shock delivery' users press the button to deliver shock.



### VOICE RECORDING

Enables a check to see what kind of on-site checks were done. When the device turns on the A16 records automatically.



### PADS QUALITY CHECK

Checks connectivity, expiry date or whether they have already been used and damaged.

- Analysis to unit charged and ready in 10 seconds or less
- Lightweight (weighs less than 2kg)
- Super-fast shock delivery time
- 3 language selection – English, Chinese (Mandarin), Korean.
- Rated to IP55
- Instant switch between adult and child mode
- Battery Capacity 3000mAh/15V (200 shocks or 6hr monitoring)

- Voice prompt, Action Icon (LED) Indicator, Status LCD
- Status indicator showing the battery level, self-test (pass or fail), temperature, condition and pads
- Self-test: POST (Power on Self-Test), Periodic Test (Daily, Weekly, Monthly), Pad Quality Check, BIST (Battery Insertion Self-Test)
- SD card, IR communication support
- Optional event review computer software
- Meets 2015 AHA / ERC guidelines

## EASY TO USE



**1** • Open the cover.  
• Turn on the unit




**2** • Select language.  
• Select adult or child.  
• Attach pads.



**3** • Unit will analyse and charge.  
• Semi-automatic - press the flashing shock button.  
• Fully-automatic - shock is delivered without the need to press the button.  
• Start CPR

### OPTIONAL EVENT REVIEW SOFTWARE

The event data is stored on the SD card. When the AED doesn't have an SD card in it, or is unreadable, corrupted or damaged, the event data will be stored in the internal memory.

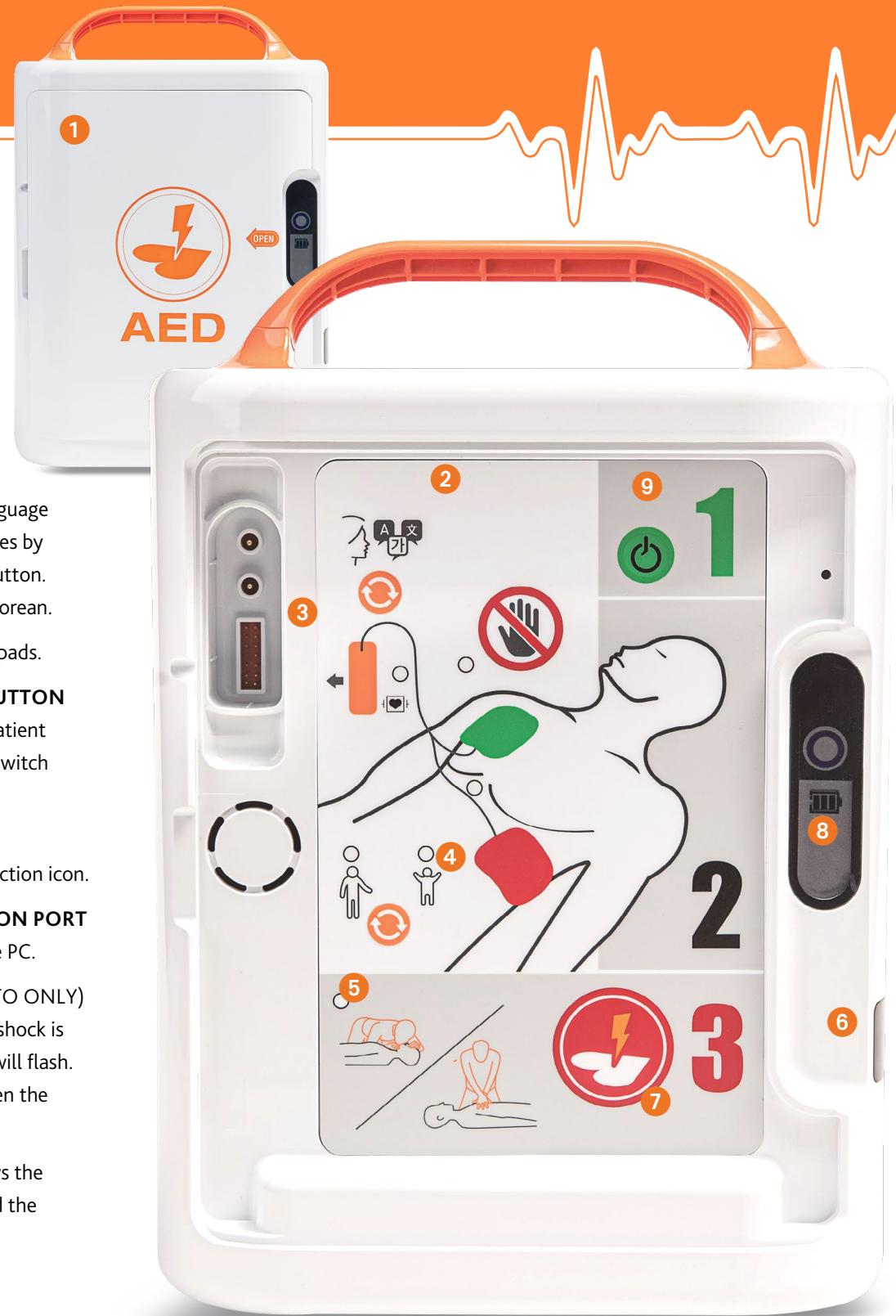

## Q and As

**Q: Can this be used for both adults and children?**

A: Yes, and the mode can be changed using the adult to paediatric mode button on the main interface.

**Q: What languages are the unit's voice prompts in?**

A: There are 3 languages loaded onto the A16, these are English, Chinese (Mandarin), Korean.

**Q: What temperature can the unit be stored at?**

A: Between 0° and 43° Celsius.

**Q: What IP is this rated to?**

A: The A16 unit is rated to IP55.

**Q: How many years of warranty does the unit have?**

A: The A16 has 5 years of warranty.

**Q: What is the difference between the two types of pads?**

A: Code 2902 are the normal adult/paediatric pads. Code 2903 include a "pads quality" function, which checks whether there is a problem with the pads, for example, expired pads, damaged pads or the recommended replacement of the pads.

**AED**

## AED SPECIFICATIONS

### DEFIBRILLATION ELECTRIC SHOCK

Waveform:	Biphasic Truncated Exponential (BTE) waveform (Impedance compensation)
Energy:	Adult: 170 to 195J (±5%) Paediatric: 44 to 51J (±5%)
Operating Mode:	Semi-Auto, Fully-Auto

### ECG

Lead:	II (RA, LL)
Patient Impedance:	25 to 200 ohm
Detection V/F	greater than or equal to 0.2mV

V/T      Adult: greater than or equal to 150 bpm

### INDICATION

#### CONTROLS

Standard:	Power button, shock button, patient mode switch button, select language button.
-----------	---

#### INDICATORS

Visible:	Action icon, Status LCD (AED, battery, temperature and pads status).
Audible:	Audio speaker (voice prompt, CPR indication). Beep (CPR indication, power on, critically low battery, self-test fail, alarm of abnormal operation).

### PHYSICAL

Dimensions:	286.5H x 200W x 90mmD
Weight:	Approx. 1.95kg (including battery, excluding pads).
Ingress Protection:	IP55

### SELF TEST

Cycle:	Every 24 hours, 1 week, 1 month. Power on self test, battery insertion self test.
Test Result:	Status LCD displays 'O' / 'X'.

### DATA BACKUP AND COMMUNICATION

Standard	Internal memory, SD card slot, IR communication port.
----------	--

## ACCESSORY SPECIFICATIONS

### PADS

#### ADULT/PAEDIATRIC

Shelf Life:	36 months from manufacture date (ref 2902)
-------------	--

Electrodes:	Disposable pads
-------------	-----------------

Placement:	Adult: Anterior-lateral Paediatric: Anterior-posterior
------------	---

Minimum active gel area:	80cm <sup>2</sup> +/- 5%
--------------------------	--------------------------

Cable Length:	Approx. 1.8m
---------------	--------------

### ENVIRONMENTAL CONDITIONS

Temperature:	Operation: 0 to 43°C (32 to 109.4°F) Storage: 0 to 43°C (32 to 109.4°F)
--------------	--

Relative Humidity:	5 to 95% RH (non-condensing)
--------------------	------------------------------

### BATTERY

Type:	LiMnO <sub>2</sub> , disposable, long-life primary cell
-------	---

Voltage/Capacity:	15V, 3000 mAh
-------------------	---------------

Shelf Life:	2 years from date of manufacture (in original packaging).
-------------	--

Standby Life:	5 years from date of manufacture (inserted in the AED)
---------------	---

Discharge:	A minimum of 200 shocks (except the CPR period between the defibrillation therapy) or more than 6 hours of operating time at 20°C
------------	---


### ENVIRONMENTAL CONDITIONS

Temperature:	Operation: 0 to 43°C (32 to 109.4°F) Storage: 10 to 25°C (50 to 77°F)
--------------	--

Relative Humidity:	5 to 95% RH (non-condensing)
--------------------	------------------------------

*The pad configuration must always match device configuration.*

*Information correct at time of going to press. Subject to change without notice.*

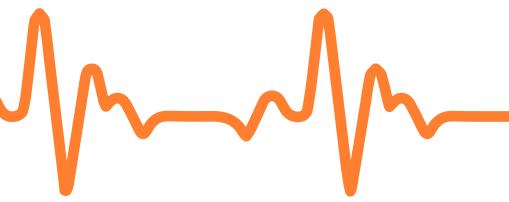

## MEDIANA AED A16



<b>A16M-DF-0E</b>	HeartOn A16 AED Auto + Multilingual Voice Prompt English, Chinese (Mandarin), Korean
CONTENTS	
1	Mediana Semi-automatic HeartOn AED A16
1	Operators Manual
1	Battery Pack
1	SD Card
1	Set Adult/Paediatric Pads
20cmH x 28.7cmW x 9cmD	
1	SIZE:
1	UNIT:
1	MOQ:



<b>A16M-DS-0E</b>	HeartOn A16 AED Semi Auto + Multilingual Voice Prompt English, Chinese (Mandarin), Korean
CONTENTS	
1	Mediana Fully-automatic HeartOn AED A16
1	Operators Manual
1	Battery Pack
1	SD Card
1	Set Adult/Paediatric Pads
20cmH x 28.7cmW x 9cmD	
1	SIZE:
1	UNIT:
1	MOQ:



<b>A16M-DF-VQ-0E</b>	HeartOn A16 AED Auto + Voice Recording + Quality Check English, Chinese (Mandarin), Korean
CONTENTS	
1	Mediana Fully-automatic HeartOn AED A16
1	Operators Manual
1	Battery Pack
1	SD Card
1	Set Adult/Paediatric Pads
20cmH x 28.7cmW x 9cmD	
1	SIZE:
1	UNIT:
1	MOQ:



<b>2902</b>	Mediana A16 HeartOn AED Adult / Paediatric Pads 16B
SIZE:	Adult/Paediatric
UNIT:	1
MOQ:	1



<b>2905</b>	Mediana A16 HeartOn AED Battery Pack
SIZE:	A16
UNIT:	1
MOQ:	1



<b>2906</b>	Mediana A16 HeartOn AED SD Card
SIZE:	A16
UNIT:	1
MOQ:	1



<b>9220</b>	Spectra A16 AED First Aid System
SIZE:	75cmHcm x 34.7cmW x 10.5cmD
UNIT:	1
MOQ:	1



<b>2910</b>	Mediana T16 HeartOn AED Trainer
CONTENTS	
1	Mediana HeartOn AED Trainer T16
1	Operators Manual
1	Remote Control
1	Set Trainer Pads
20cmH x 28.7cmW x 9cmD	
1	SIZE:
1	UNIT:
1	MOQ:



<b>2911</b>	Mediana T16 HeartOn AED Trainer Pads
SIZE:	
UNIT:	1
MOQ:	1

## AED CONSUMABLES



<b>A2780</b>	AED Prep Kit in Extra Small Aura III Box
SIZE:	20.5cmH x 20cmW x 6.5cmD
UNIT:	1
MOQ:	1



<b>2877</b>	AED Prep Kit in Black Soft Pouch
SIZE:	15.6cmH x 10.8cmW x 5.4cmD
UNIT:	1
MOQ:	1



<b>2102</b>	AED Alarmed Metal Storage Cabinet White/Glass/Alarmed
SIZE:	Internal: 43.4cmH x 43.4cmW x 21.7cmD External: 44.4cmH x 44.4cmW x 22.7cmD
UNIT:	1
MOQ:	1



<b>2879</b>	Universal AED Bracket Green
SIZE:	12.3cmH x 19cmW x 9.5cmD
UNIT:	1
MOQ:	1



<b>3098</b>	AED Alarmed First Aid Cabinet White/Glass/Lockable
SIZE:	External: 46cmH x 30cmW x 14cmD
UNIT:	1
MOQ:	1



<b>7900</b>	LifeVault AED Cabinet Locked
SIZE:	xxxx
UNIT:	1
MOQ:	1



<b>7901</b>	LifeVault AED Cabinet Unlocked
SIZE:	xxxx
UNIT:	1
MOQ:	1



# MEDIANA HEARTON AED A16

Automated External Defibrillator

[medianadefib.co.nz](http://medianadefib.co.nz)

New Zealand

Tel : +64 22 623 5237

email: [kiaora@reliancemedical.co.nz](mailto:kiaora@reliancemedical.co.nz)

Website: [www.reliancemedical.co.nz](http://www.reliancemedical.co.nz)

Australia

Tel : 1300 806 120

Email: [rmau@reliancemedical.au](mailto:rmau@reliancemedical.au)

Website: [www.reliancemedical.co.nz](http://www.reliancemedical.co.nz)



Mediana HeartOn Defibrillators are distributed  
by Reliance Medical Limited.

